(climsom Zen)

EYE MASSAGER

Mode d'emploi

User manual

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(Find a quiet place in the noisy world)

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Bonne(s) séance(s) !

Thank you for choosing Climsom Zen Eye Massager.

Bring convenience and comfort to your life thanks to this exceptional way to relax, take care of your eye contour, improve your sleep, relieve eye fatigue and prevent headaches.

Before operating this unit, please read these instructions completely to enjoy your massage and take full advantage of its excellent performance.

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PRECAUTIONS

Notes

Make sure to carefully read and follow the operation guidelines of this user manual before using the product. The company cannot be held liable for any consequence arising from incorrect usage or failure to abide by this user manual.

Don't use the product after wearing make-up as the product may mess your make-up and cosmetics may pollute the product surface. In case the product is polluted, wipe it with clean water.



Make sure you are not wearing contact lens or glasses and remove accessories that may cause discomfort during use.

Warning

- Patients who had eye surgery and heart disease should not use the product.
- Children and pregnant women should not use the product.
- Don't use the product if you are suffering from stroke, high blood pressure or low blood pressure.
- Don't use the product if you have eye symptoms such as glaucoma, cataract, amotio retinae, pinkeye.
- Don't use the eye massager during driving.

- Don't use the eye massager at places subject to very high temperature or strong humidity.
- Avoid using the eye massager after drinking alcohol.
- Users subject to sensitive skins, facial allergic dermatitis, facial cosmetic surgery are advised to use the product under the guidance of doctors.
- The product can be heated. Users who are not sensitive to heat should be careful when using it.
- The product is a personal care massage product and cannot replace medicines and medical devices.
- As the product contains lithium battery, don't place it in high temperature environment.
- Don't massage parts other than your eyes area with the product.
- Don't disassemble the product without authorization.
- Immediately stop using the product in case of discomfort during usage.

PRODUCT DISCOVERY

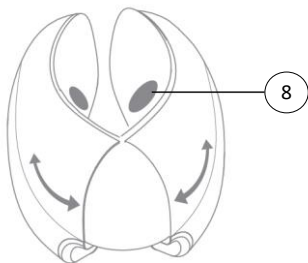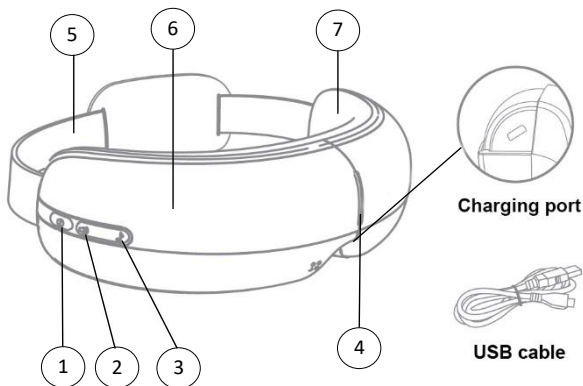1. On/Off button
2. Volume button
3. Music button
4. Light indicator
5. Elastic band
6. Upper cover
7. Leather cover
8. Air bag

Overview of Eye Massager various functions:

- **5 message modes** with mix of air pressure, hot compress, vibration.
 - > *Long press On/Off button to turn on or turn off the massager*
 - > *Short press On/Off button to switch modes*
- **3 sound atmospheres:** seaside, nature and Zen music. You can also turn the music off.
 - > *Short press Volume button to adjust the volume*
 - > *Short press Music button to switch sound atmosphere*
 - > *Long press Music button to turn off music*
- **Bluetooth mode** to enjoy your own music
 - > *Default enabled Bluetooth when Eye Massager is switched on*
- **Adjustable elastic band** for a more or less intense massage
- **Charging battery** to use your wireless Eye Massager (USB cable included)
 - > *Battery indicator above Music button (can be seen when device is turned on): flashes when device needs to be charged*
- **15 minutes default working time**
 - > *Long press On/Off button to restart Eye Massager*
- **Easy transport** (the massager can be folded) to bring you comfort wherever you go

USE

Preparation

Your massage device works with USB rechargeable battery. If the battery is discharged, you will not be able to switch on your massager. **Please charge the product before use** (see page 23 "Charging") by connecting your device to a USB port (computer, USB wall socket, AC adapter...) thanks to the USB cable provided. During charging, the battery indicator, above Music button, is red. Charging time is between 3 and 5 hours. When the light turns green, the massager is charged and ready for use.

How to use

1. Put the device over the eyes. Adjust the headband. You can use the device sitting or lying down.
2. **Press On/Off button for two seconds to turn on the massager** with the voice prompts "*Hi*" then "*Mode one*" which is default enabled message mode.
3. Short press On/Off button to switch message mode (see next page "Message mode selection").
4. After 15 minutes the message will turn off with the voice prompting "*Thanks, bye!*".

To extend usage time, reboot the massager. To shorten your session, press On/Off for 2 seconds to turn the unit off.

Massage mode selection - The product has five modes:

	 Air pressure	 Hot compress	 Vibration
Mode 1	x	x	x
Mode 2	x	x	
Mode 3	x		x
Mode 4		x	
Mode 5	x		

What do you need?

Active mode: select Mode One.

Soothing mode: select Mode Two.

Eye care mode: select Mode Three.

Sleep mode: select Mode Four.

Relaxation mode: select Mode Five.

Don't hesitate to test every mode to enjoy the best one regarding your need and preference.

You can see massage symbols of mode in progress above Volume and Music buttons when the unit is switched on.

For example, if you select Mode 1, you will see the following symbols: 

Volume control - The default volume for booting is level 3. Short press Volume button to switch to 3-4-5-off-1-2-3... and so forth.

Music selection - You can choose between 3 musics: seaside, nature or Zen. By default, seaside music is on when booting followed by nature music and Zen music. Short press Music button to switch (in the order seaside/nature/Zen) and listen only one music during whole massage. Long press to pause the music and enter silent mode. To listen your own music, use Bluetooth.

Bluetooth - For first time use, please connect Eye massager with the audio equipment of your choice (smartphone, computer or any other device containing your music and having a Bluetooth connection):

1. Turn your Eye Massager on (long press On/Off button). Bluetooth is enabled by default.
2. Activate Bluetooth on your audio equipment.
3. Connect to the Eye Massager Bluetooth named "E3 Eye Massager".

4. Once connected, the voice prompts "*Connected*": you can play your music.

If you disconnect Bluetooth from your audio equipment, the voice prompts "*Disconnected*" and your music stops playing. Thereafter, Bluetooth connection will be automatic when the device is switched on and Bluetooth activated on your smartphone or other media.

Charging

- Before first use, please charge the product using the USB cable (provided) to connect your device to a computer or using an AC adapter (not provided).
- When the battery is low, the massager will prompt "*Battery low*" and battery indicator above Music button will flash. Please fold the product to expose charging interface and insert the USB cable provided. Plug it into a USB port (from a computer for example). The charging time of the built-in battery ranges from 3 to 5 hours (depending on USB port - 5V2A is faster than 5V1A).
- During charging, battery indicator (located above Music button) is red and the product cannot be used (the product initiates self-protection procedure to stop working automatically when it's charged).
- After a full charge, battery indicator becomes green: you can use your device again. As an indication, a full battery allows 6 to 8 sessions.

After the product is fully charged, make sure to remove the power adapter to avoid long-time charging which reduces product life.

- During usage, if the voice prompts low battery alarm (“*Battery low*”), we recommend to stop your session and start it again after a full charge so as to preserve battery life.

Daily maintenance

- Use clean water and wet cloth to wipe the dirty surface. Don’t use alcohol and other volatile liquid to wipe. The product cannot be disassembled and put in water for washing.
- During cleaning process, prevent water or others liquids from entering the machine which could damage it and create a short circuit.
- For long-term storage, cover the product to avoid dust: you can put in its box to protect it.

Troubleshooting

Unable to boot

- Low power: fully charge the massager.
- Charging in progress: use after charging.

Unable to charge or no charging indication

- The USB cable, massager and power adapter are not well-

connected: re-connect after unplugging. Try to connect to another USB port.

- The USB cable is damaged: change the USB cable.

WARRANTY POLICY

During the warranty period, you can enjoy the services of replacement and refunding according to the following policy:

1. 14 days withdrawal period from the day after receipt of the goods in case of distance selling and depending on the local regulations. Return costs for problems other than quality are at the customer's expense. Please return the product in its original packaging with all accessories. Please contact your seller before returning the product.
2. Climson guarantees the purchaser against any manufacturing fault for 2 years from the date of purchase. The implementation of the guarantee will result the only replacement of the defective product to the exclusion of any other repair.

What is not covered by the warranty:

1. The product has been tampered with, repaired or modified by non-authorized personnel;
2. Damage (accidental or otherwise) to the product that does not impact the product's operation, such as without limitation to rust, change in color, texture or finish, wear and tear, and gradual deterioration;

3. The damage of products caused by accident, abuse, natural disaster, intentional or accidental misuse, fire, neglect or improper maintenance and use under abnormal conditions.

SPECIFICATIONS

Model: Climsom Zen Eye Massager – E3

Battery capacity: 1100mAh

Rated input power/voltage: 5W/5V

Input interface: Micro USB

Charging time: 3-5 hours

Size: 8.27 inches (6.30 inches once folded)

Weight: 0.78 pound

Material: plastic ABS + sea yarn texture fabric for upper cover, skin-friendly protein leather for leather cover

Included: Eye Massager + USB cable + user manual + storage box

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*Discover CLIMSOM other products for your well-being,
beauty and health on **www.climsom.com***

We will be glad to help you:
+33 (0)2 85 52 44 74
contact@climsom.com

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