

PRECAUTIONS

- This product is designed for your comfort and for a domestic use. If there is any doubt about the heat or coolness use in your situation, consult a doctor.
- In cases of sickness or discomfort, stop using CLIMSOM INTENSE and consult a doctor.

ADVICE

Use the intense freshness to relieve heavy legs and relieve restless legs.

Use the intense heat on low back pain, lumbago and other muscle contractures.

The coolness or heat will be felt even through clothing.

For even more freshness or warmth felt, use the CLIMSOM INTENSE mattress topper without its cover.

Put yourself on the "soft" fabric surface of the CLIMSOM INTENSE mattress when you want a more moderate temperature.

You can program a delayed start or an early stop. Refer to your thermo-controller instructions manual.

COMPOSITION

Cover – plastic surface	PVC / NYLON
Cover – fabric surface	80% POLYESTER – 20% COTON
Inside mattress	100% PVC
Connection pipe & hose	100 SILICONE

Discover our other wellness, bedding products and in particular the Climsom

Canicule, large version of Climsom Intense,
on our website www.climsom.com.

Contact us by email at contact@climsom.com
or by phone on +33 (0)2 85 52 44 74.

(climsom) *intense*

Thermotherapy Mattress Topper




To work, the CLIMSOM INTENSE mattress topper must be connected to a CLIMSOM thermo-controller.

The CLIMSOM INTENSE mattress topper achieves the coolest or warmest temperatures, because of its PVC material which improves thermal transfer and because of the large amount of water flowing through the PVC “grid”.


Read carefully this manual before use



DESCRIPTION

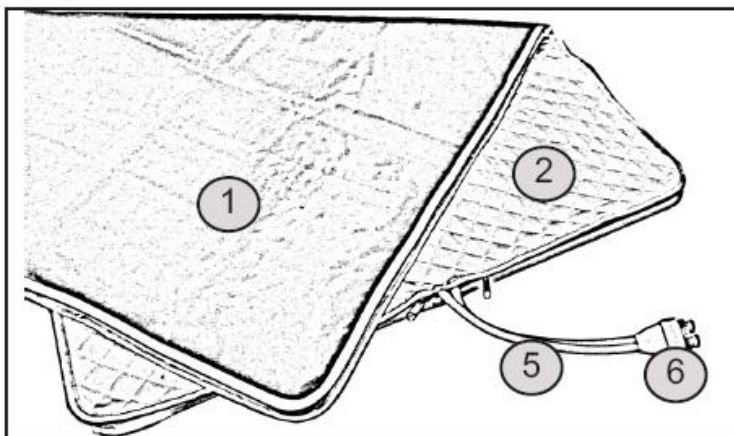

- 1 – « Intense » thermo-controlled surface
- 2 – « Soft » thermo-controlled fabric surface
- 5 – Connection pipe
- 6 – Connector

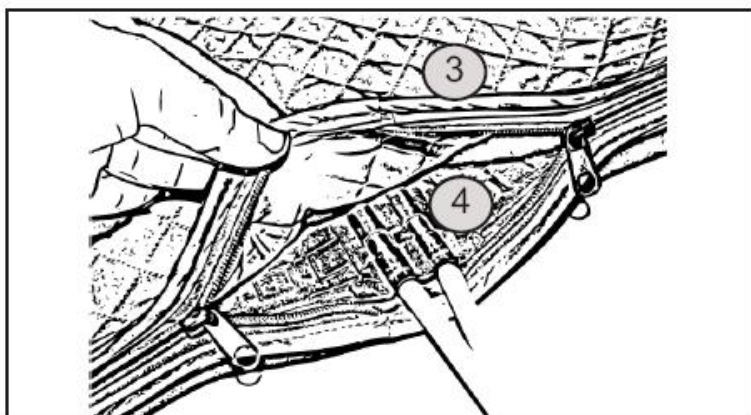

- 3 – Machine cashable cover
- 4 – Inside PVC mattress topper

INSTALLATION

1 – Place your CLIMSOM INTENSE mattress topper between your mattress and your fitted sheet, or directly on your fitted sheet. Place the "intense" plastic thermo-controlled surface (1) face-up for direct contact with the user or choose the "soft" fabric thermo-controlled surface (2) if you want less intense heat. You can also remove the mattress topper cover for direct contact with the inner PVC mattress topper (4).

2 – Connect your CLIMSOM INTENSE mattress topper to your CLIMSOM thermo-controller with the 1,20m hose included. Be sure each connector is correctly plugged in. You must hear a « click ».

3 – Fill the water tank of your CLIMSOM thermo-controller, ideally with demineralized or distilled water, or with low mineral content and low lime.

4 – Turn on your CLIMSOM thermo-controller and set the temperature on the screen or with the remote. To start, try moderate temperatures and increase or decrease gradually to find your personal comfort temperature.

The CLIMSOM INTENSE makes it possible to reach temperatures between 15 ° C and 48 ° C, these temperatures are given as an indication as they depend on many factors (volume of the room, surrounding temperature, sleeper temperature, room ventilation, sleeper morphology, thermal inertia of the mattress...).

MAINTENANCE

The silicone pad (4) can be cleaned with soft soap and sponge, without any detergent. The cover (3) can be machine washed (gentle cycle, 30°C without bleach water), after removing the inside PVC pad (4).

- Let dry in open air. Do not use tumble dryer.
- Do not expose to sunlight.
- Do not iron.
- Do not dry clean.
- Do not mark the folds when you store your CLIMSOM INTENSE.
- For storage, empty water from the pad by pressing on the "drop-stop" springs connector thanks to the drain key included with your CLIMSOM thermo-controller, and place the product in a dry place.
- Due to temperature spread, there might be condensation: air out your bed.

